

# INTRODUCING

World's first dual function plasma freezing solution

>> INNOVATIVE >> FAST >> AFFORDABLE  
>> MADE FOR INDIA.

# CRYOPRES DUO STORE




REMI's **Cryopres Duo Store** combines a **Contact Shock Freezer** and a **-40°C Deep Freezer** in one vertical unit, **revolutionising Plasma Freezing & storage.**

Designed specifically for blood banks, this innovation ensures that **plasma is frozen within one hour preserving essential coagulation factors and meeting pharmaceutical grade plasma quality requirements.**

# CRYOPRES DUO STORE

**A Contact Shock Freezer for Every Blood Bank**

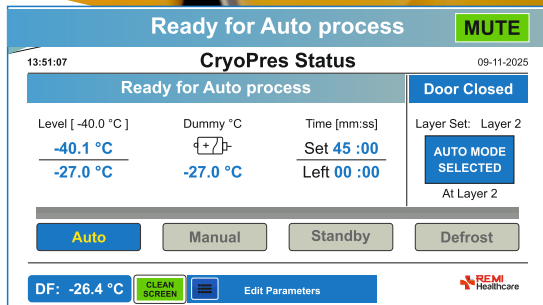

### FREEZING PLATES

#### CryoPres Duo Store

Experience the unique dual function design which is powerful combination of contact shock freezer & -40 deg C Deep Freezer in one compact unit

#### Smart, Intuitive Operation

- Touch Screen Interface for effortless control of both freezing and storage functions.
- Real-time display of date, time, temperature, and operating status.
- Independent & alternate screen display for Contact Shock Freezer and Deep Freezer.
- Multi-alarm system for High/Low Temperature, Door Open, Power Failure, and Freezing in Progress alerts.
- Hot gas based Auto defrost.



#### Precision + Safety

- The Contact Shock Freezer door auto-locks during freezing cycle to ensure consistent thermal performance.
- Freezing Capacity: 12 Plasma Bags frozen to the core in just 45-50 minutes.
- Storage Capacity: 336 Litres, ideal for long-term plasma preservation.
- 7-day circular chart recorder for reliable temperature documentation of the storage compartment.

#### Built for Performance

- Inner & Outer Body: High-grade Stainless-Steel construction for durability and hygiene.
- Dual Function, Dual Assurance: Even when the Contact Shock Freezer is not in use, the display seamlessly switches to show the storage temperature.

# CRYOPRES

## DUO STORE

### A Contact Shock Freezer for Every Blood Bank